

M7

Hand-carried Color Doppler Diagnostic Ultrasound System

Musculoskeletal Solutions


- iZoom: One key magnification, instantly expands the image to full screen for better visualization in distance
- B-steer: Steered beam line for difficult inspection area
- Trapezoid Imaging: Extends field of view for more anatomy information
- iScape View: Real-time panoramic imaging, enlarges view for large anatomy structures
- iNeedle: Needle enhancement visualization for better observation of the needle and tip
- Dedicated needle guide brackets for interventional procedures
- Abundant transducers for different scenarios




Needle Guide Bracket




Transducers

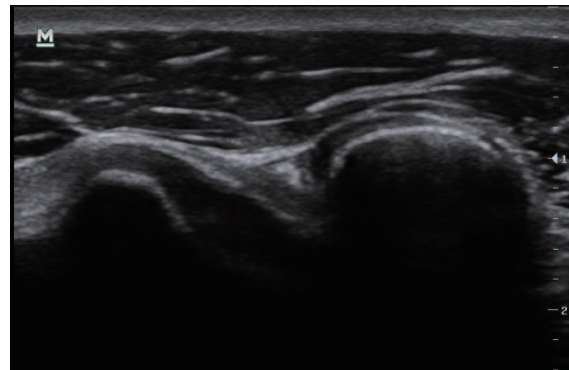


Acromioclavicular Joint

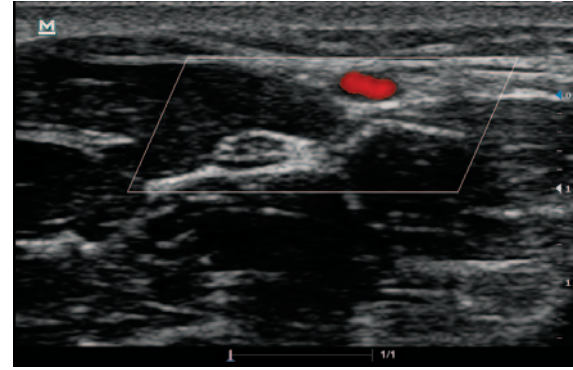

Median Nerve & Radial Artery

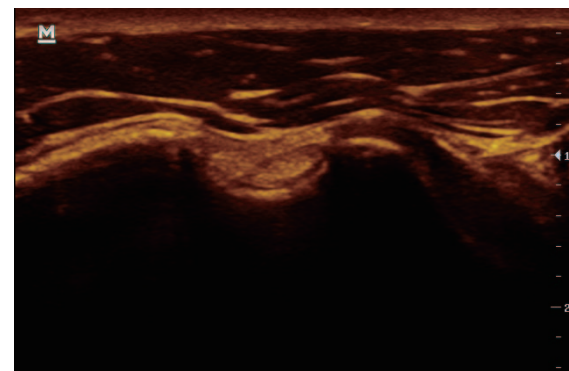

Right Shoulder

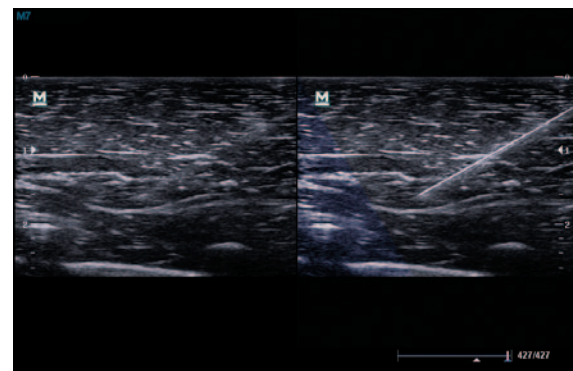

iNeedle



M7

Hand-carried Color Doppler Diagnostic Ultrasound System

Expert in Musculoskeletal Imaging

Excellent Image Quality

- 3T Transducer Technology: provides better spatial resolution and noise reduction
- iBeam: Spatial compounding imaging, enhances boundary display
- iClear: Speckle suppression imaging, provides better contrast resolution for tissue differentiation
- Phase Shift THI: provides better contrast resolution
- High frequency linear transducer: broadband and multi-frequency for high resolution imaging of muscles, tendons and ligaments

Mobility and Ease of Use

- Short boot-up time with batteries for point-of-care
- Anti-shock magnesium case
- Mobile trolley with extra monitor: expands your work platform to everywhere
- Comfortable backpack and artistic travelling case
- Transducer extender: Three active transducer connectors
- iStation: Patient information management system for data collection and storage
- iTouch: Optimizes the parameters of image by single push

mindray

Mindray is listed on the NYSE under the symbol "MR"
Mindray Building, Keji 12th Road South, High-tech Industrial Park,
Nanshan, Shenzhen 518057, P.R. China
Tel: +86 755 8188 8998 Fax: +86 755 26582680
E-mail: intl-market@mindray.com Website: www.mindray.com



MINDRAY is a trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. Specifications subject to changes without prior notice.
© 2012 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved.
P/N: ENG-M7-MI-210285X2-20120317

mindray
healthcare within reach